



SAB CLASSIC

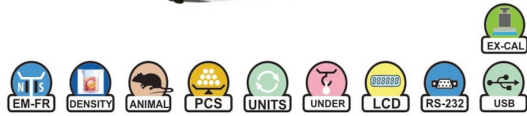

ANALYTICAL & PRECISION BALANCE

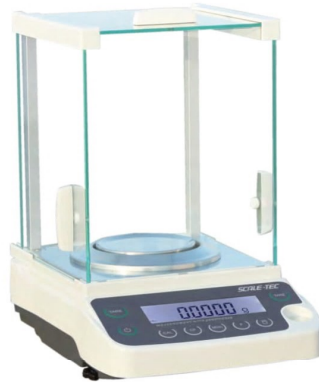

Get extensive high resolution even when you weigh high capacity with latest Monolithic Block Technology.

MODEL	CAPACITY (G)	RESOLUTION (G)	PAN (MM)	REPEATABILITY (G)	LINEARITY (G)	RESPONSE TIME (SEC)
External Calibration						
SAB 124CL	120	0.0001	90	± 0.0001	± 0.0002	< 3
SAB 224CL	220	0.0001	90	± 0.0001	± 0.0002	< 3
SAB 223CL	220	0.001	110	± 0.001	± 0.0015	< 3
SAB 323CL	320	0.001	110	± 0.001	± 0.0015	< 3
SAB 423CL	420	0.001	110	± 0.001	± 0.0015	< 3
SAB 523CL	520	0.001	110	± 0.002	± 0.002	< 3
SAB 623CL	620	0.001	110	± 0.002	± 0.002	< 3

Option Accessories



Features

The New Price & Best in Performance Class Analytical & Precision Balance.

- * Single Touch Key External calibration, ensures the highest accuracy at the press of a key.
- * Support: Weighing, Counting, Percentage, Animal Weighing. Optional: Density measurement
- * 18 standard units, with all-character mark.
- * Strong anti-shaking monolithic electromagnetic for highest accuracy.
- * Standard 1500V isolation USB port, ease to connect with PC.
- * User selectable filter, stability indicator and auto-zero tracking for use in difficult environment.
- * Large, sturdy and nearly frameless, all-glass draft shield chamber for optimal viewing of weighing chamber.
- * The draft shield doors and all stainless steel parts of the weighing chamber base are easy to remove in just a few steps for cleaning.
- * Backlit, high-contrast display with 15 mm digits; exceptionally easy to read under any room lighting conditions.

SAB CLASSIC INCAL



BASIC ANALYTICAL BALANCE - INCAL



Get extensive high resolution even when you weigh high capacity with latest Monolithic Block Technology.

MODEL	CAPACITY (G)	RESOLUTION (G)	PAN (MM)	REPEATABILITY (G)	LINEARITY (G)	RESPONSE TIME (SEC)
Internal Calibration						
SAB 124CL IN CAL	120	0.0001	90	± 0.0001	± 0.0002	< 5
SAB 224CL IN CAL	220	0.0001	90	± 0.0001	± 0.0002	< 5

Option Accessories



Features

The New Price & Best in Performance Class Analytical Balance.

- * Built-in motorized calibration weight, ensures the highest accuracy at the press of a key.
- * Support: Weighing, Counting, Percentage, Animal Weighing. Optional: Density measurement
- * 18 standard units, with all-character mark.
- * Strong anti-shaking monolithic electromagnetic for highest accuracy.
- * Standard 1500V isolation USB port, ease to connect with PC.
- * User selectable filter, stability indicator and auto-zero tracking for use in difficult environment.
- * Large, sturdy and nearly frameless, all-glass draft shield chamber for optimal viewing of weighing chamber.
- * Backlit, high-contrast display with 15 mm digits; exceptionally easy to read under any room lighting conditions.