


CARBON & SULPHUR APPARATUS (STROHLEIN APPARATUS)



- The furnace is designed to accommodate combustion tubes of dimensions 19mm x 24mm x 600mm. Tubular combustion furnace unit, heated by 2 Nos. of Silicon Carbide heating elements suitable for 1300 C.
- Control panel with Ammeter, Voltmeter, Temperature is controlled by Thyristor Based Digital Temperature Controller Cum Indicator **from 50 deg C to 1300 deg C. +/-1 deg C. with Cr/Al Thermocouple.**
- Glass Parts of borosilicate Glass are mounted on an attractive Fibre stand.
- Supplied Complete Kanthal Scoop, Clamps for heating Elements, Rubber tubing, Combustion Boats, Combustion Tube and Silicon Carbide Heating Elements.

www.everflowscientific.com

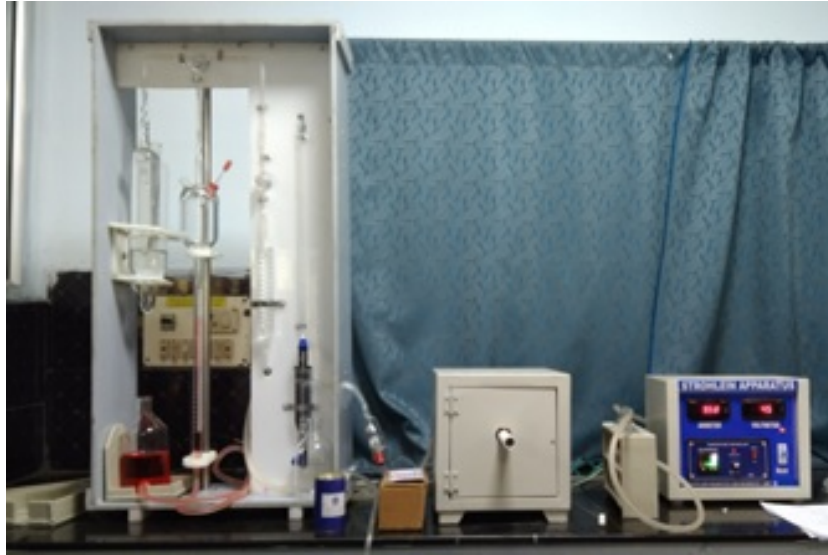

EVERFLOW SCIENTIFIC INSTRUMENTS

No. 231 / 298, Kilpauk Garden Road, Kilpauk, Chennai - 600 010.

Mob: 9444259624, Ph : 044-26646066, 26645708 e-mail: everflowin@gmail.com

GSTIN: 33AFBPG5986P1ZX